

# veinlite PEDI

## For Children, Babies and Neonates

### No More Guesswork - No More Probing



### Proven Veinlite® Technology

Veinlite PEDI is our smallest and most compact version of the clinically proven and very successful Veinlite LED device which was reviewed as *"an invaluable new addition to our kits."* (Jeffery Lindsey - JEMS 2005)

### Invaluable for Vein Access

Take the stress out of finding veins in your smallest patients! The Veinlite PEDI is optimized to find veins using the traditional through-the-body transillumination technique and also by the patented\* side-transillumination method which can be used anywhere on the body, even on the scalp. With the Veinlite PEDI, veins can be mapped to avoid valves or bifurcations and find the best place for safe access. The ingenious C-shaped design helps secure the vein, provides traction on the skin, facilitating easy needle insertion.


### Snapshot

- ♦ LED colors optimized for vein imaging
  - Orange LEDs (8) for superficial veins\*\*
  - Red LEDs (4) for deep veins, dark skin
- ♦ Disposable clear plastic cover
- ♦ Small 3 volts Lithium Battery
- ♦ Weight 40 gm (1.5 oz) with battery
- ♦ Size 6.7 x 4.2 x 2.4 cm
- ♦ One year warranty



### Unique C-Shaped Design

With Veinlite PEDI's C-shaped design and 12 light emitting diodes (LEDs) in two colors, superficial veins can be seen and accessed anywhere on the body. No more guess work and no more probing.

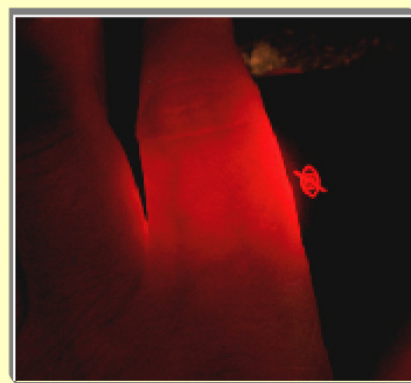

### Snap-on Neonate Adapter

An innovative adapter\*\* allows six LEDs to shine through a small opening for through-the-body transillumination of limbs.

## TransLite

medical transillumination

8410 Hwy 90A, #150, Sugar Land, TX 77478  
 www.veinlite.com info@veinlite.com  
 Tel: 281 240 3111 Fax: 281 240 3122

\* - US Patent Number 5,146,923, \*\* - Patent Pending, © Registered Trademark of TransLite LLC

VPED0806